



Gentle for the newborn,
efficient for you.

DRÄGER
JAUNDICE METER
JM-105

Assessing jaundice risk – without the trauma



While neonatal jaundice – or hyperbilirubinemia – occurs in almost all babies, it is a cause of concern for caregivers and anxiety for families. If not treated in time, neonatal jaundice can lead to permanent brain damage¹.

Screening for jaundice by visual assessment can lead to an overestimation of risk, which means unnecessary lab tests. This blood draw (total serum bilirubin testing (TsB)) requires a heel prick, which is painful for the baby and costly for the hospital. Visual assessment can also lead to an underestimation of risk, which could result in a failure to obtain necessary lab tests.

That's why transcutaneous bilirubin testing (TcB) has grown rapidly in the past few years as a standard practice in hospitals to identify at-risk infants. But even TcB has presented challenges in terms of human error and time consumption. Until now.

¹ Bhutani VK, Johnson LK, Keren R. Diagnosis and management of hyperbilirubinemia in the term neonate: for a safer first week. *Pediatr Clin North Am* 2004;51(4):843-61.



NEXT-GENERATION JAUNDICE SCREENING FROM DRÄGER

Building on more than 30 years of proven technology², the Dräger Jaundice Meter JM-105 is a non-invasive transcutaneous bilirubinometer that measures the yellowness of subcutaneous tissue in newborn infants as young as 24 weeks gestational age³ who have not undergone a transfusion or phototherapy. It provides a visual digital measurement that correlates with serum bilirubin results. Effective on heterogeneous populations, the JM-105 can be used in both hospitals and outpatient settings.

The JM-105 improves the process of jaundice screening in several ways. The risk of infection is reduced since screening with the JM-105 is non-invasive and drawing blood is not necessary. The device also reduces risk of human error because it can scan nurse and baby ID information with a bar code scanner and can automatically measure, save and transfer data. Not only does this eliminate manual transcription – which can be error prone and time consuming – but it also results in quicker access to information for decision making.

Screening with the JM-105 is cost effective because it reduces the frequency of lab tests. And because the device has a reusable probe, expensive disposables are not needed.

Because bilirubin levels peak two to four days after birth – when most babies have gone home – jaundiced babies have to be readmitted for treatment. Since the JM-105 can detect this risk before babies leave the hospital, jaundice can be treated immediately, which can reduce readmissions and lengths of stay⁴.

² Transcutaneous Bilirubinometry; Daniele De Luca, William Engle & Greg Jackson, Published by Nova Science Publishers, Inc. New York. 2013

³ Dr. Jackson & Dr. Engle, Evaluation of transcutaneous bilirubinometry in preterm neonates, 2009, Journal of Perinatology 29

⁴ Peterson JR, Okorodudu AO, Mohammad AA, Fernando A, Shattuck KE. Association of Transcutaneous Bilirubin Testing in Hospital with Decreased Readmission Rate for Hyperbilirubinemia. Clinical Chemistry 2005;51

Introducing a **better** way to deliver care to your newborns

Patient flagging

Allows you to easily record and identify babies that need special attention

Lightweight, ergonomic design

Fits easily in your hand or pocket for easy transport within the hospital or between clinical offices; smooth shape is easy to clean

Integrated light checker

Light checker in docking station helps ensure optimal performance of light source

Reusable probe

Just wipe with alcohol and reuse. There are no disposable tips to fumble with, keep track of, or reorder

SIMPLE JAUNDICE SCREENING PROCESS



D-94081-2013

1. Nurse enters ID for nurse and baby with barcode scanner or touch-screen alphanumeric keypad



D-104476-2013

2. Nurse cleans probe and takes TcB measurements; no need to recalibrate for each patient



D-104476-2013

3. Nurse flags patients who may need immediate attention (!)



D-104476-2013

Bar code scanner
Provides fast, accurate entry of nurse and patient identification information

Touch screen display
Makes it easy to navigate through menu settings; choose to see results in mg/dL or $\mu\text{mol/L}$; adjust contrast and sensitivity of display


Faster measurement
Saves time because the tip is reusable and calibration is not required for every patient

Internal memory
Holds up to 100 patient files

Data transfer via HL7
Enables connectivity with electronic medical record



D-104474-2013

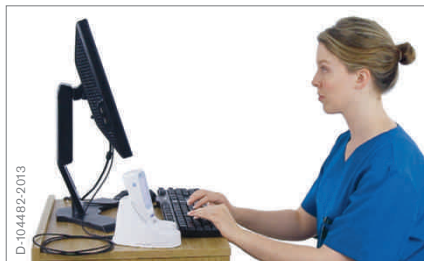

D-104482-2013



MT-9251-2003

4. Nurse docks JM-105 into docking station, which is connected to the hospital network

5. Nurse can automatically transfer data for multiple patients to your hospital's data management system to store measurement results centrally

6. Doctor reviews TcB results to decide on further medical care

Benefits patients, parents and caregivers



The JM-105 gives you consistent quality, cost-effectively delivered over the lifetime of the device. As a result you can optimize the efficiency of your jaundice program, save time and money, and deliver an exceptional standard of care for your newborns.

EFFECTIVE

- Accurately identifies at-risk infants and eliminates guesswork
- Screens infants as young as 24 weeks gestational age
- Helps you establish effective protocols and comply with the latest clinical guidelines and recommendations



EFFICIENT

- Reduces manual documentation via electronic data transfer
- Provides immediate access to information for quicker, informed decision making
- Reduces the cost of jaundice screening over lifetime of device by eliminating disposables
- Helps decrease unnecessary heel pricks, infection risk, and parent anxiety

EASY TO USE

- Makes screening fast and easy with an intuitive user interface
- Provides faster measurement since no calibration is needed for each new patient
- Offers bar code reader to enter nurse and patient ID
- Gets your team up and running quickly with supportive training tools

Simply put, the Dräger Jaundice Meter JM-105 is gentle for the newborn and efficient for you.



Neonatal competence



D-22387-2010

NEONATAL CARE

Dräger offers a comprehensive range of products to meet the complex needs of infants, families and caregivers – from innovative neonatal ventilation systems and specialized incubators to mobile monitors and state-of-the-art IT systems. Additionally, Dräger design and planning service help NICUs create an ideal ergonomic workplace where caregivers can work efficiently. Parents can nurture, families can bond, and infants can thrive.

Education and training



D-94388-2013

TRAINING RESOURCES

As part of our ongoing commitment to effective training, Dräger offers a comprehensive training video and training manual that deliver clear, concise instruction that enables you to become proficient in routine use of the Dräger Jaundice Meter JM-105 quickly. These resources allow you to train at your convenience, while helping you build confidence in using your equipment.

Service



D-32436-2011

DRÄGERSERVICE

At Dräger, we believe that quality devices deserve quality service. It is good to know that you have a partner at your side you can rely on for a wide range of service needs. We know your devices and learn your processes. We use the same quality standards for service that apply to our development of devices. As a result, your biomedical engineering department receives first class support.

BabyFirst™



D-64244-2013

BABYFIRST

BabyFirst is a preeminent online neonatal community for clinicians and parents of premature babies. With content populated by clinicians and renowned experts, this website offers families of premature babies a trusted resource for gaining a better understanding of what they can expect in the Neonatal Intensive Care Unit – with insight into common terms, procedures, equipment, post-hospital care, and more. Visit www.babyfirst.com.

Not all products, features, or services are for sale in all countries.
Mentioned Trademarks are only registered in certain countries and not necessarily in the country in which this material is released. Go to www.draeger.com/trademarks to find the current status.

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

www.draeger.com

Manufacturer:

Draeger Medical Systems, Inc.
3135 Quarry Road
Telford, PA 18969-1042, USA

Locate your Regional Sales
Representative at:
www.draeger.com/contact



REGION DACH

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION EUROPE

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION MIDDLE EAST, AFRICA

Drägerwerk AG & Co. KGaA
Branch Office
P.O. Box 505108
Dubai, United Arab Emirates
Tel +971 4 4294 600
Fax +971 4 4294 699
contactuae@draeger.com

REGION ASIA / PACIFIC

Draeger Singapore Pte. Ltd.
25 International Business Park
#04-20/21 German Centre
Singapore 609916
Tel +65 6308 9400
Fax +65 6308 9401
asia.pacific@draeger.com

REGION NORTH AMERICA

Draeger, Inc.
3135 Quarry Road
Telford, PA 18969-1042, USA
Tel +1 800 4DRAGER
(+1 800 437 2437)
Fax +1 215 723 5935
info.usa@draeger.com

REGION CENTRAL AND SOUTH AMERICA

Dräger Panama Comercial
S. de R.L.
Complejo Business Park,
V tower, 10th floor
Panama City
Tel +507 377 9100
Fax +507 377 9130
contactcsa@draeger.com