



Your customized
anesthesia workstation

DRÄGER
FABIUS® GS PREMIUM

One of the most challenging places on earth



MT-1608092009

The modern hospital environment is a challenging one. As healthcare providers continue to strive for the best possible care, the realities of increasing demands, limited budgets and rising complexity can often seem intimidating – and sometimes almost insurmountable. So when our engineers set out to design the new Fabius GS premium, they had a very simple philosophy in mind: **As practical as possible, as flexible as we can make it.**



Is there a quality anesthesia workstation out there that's flexible yet simple to use? A machine we can really rely on? One that has all the features we need and is still affordable? And one that we can customize to suit our needs?

Take a closer look at the Dräger Fabius GS premium

The Dräger Fabius GS premium now combines solid design and modular architecture with a wide range of ventilation capabilities that cover most operating room applications, delivering exceptional performance, value and flexibility. Combining time-tested design concepts with the comfort and utility of modern digital technology, it's an anesthesia workstation that is simple to use, highly efficient and ready for the future. An array of options lets you expand the capabilities of the Fabius GS premium even further to create exactly the workstation you need, letting you make the most of your investment.

PERFORMANCE FEATURES

- High contrast color monitor with Dräger user interface for easy and intuitive operation
- Supports all major ventilation modes
- New, highly maneuverable trolley with a central brake
- Integrated LED workplace illumination
- Solid, spacious design with large drawers
- Industry standard rails for additional Dräger and third party equipment
- CLIC adapter for use with Dräger sorb Clie compatible

Made with a passion for quality

HIGH PERFORMANCE VENTILATION

At its core, the Fabius GS premium features our powerful and accurate E-vent piston ventilator, a technology which requires no drive gas and delivers ICU-like ventilation performance. This means the Fabius GS premium is ready for a broad range of ventilation applications and acuity levels. All major ventilation modes can be supported, including volume control, pressure control, pressure support (PS) and SIMV/PS, providing smooth transitions from controlled modes to pressure support and spontaneous breathing. Advanced features such as dynamic compliance compensation, fresh gas decoupling and backup ventilation for PS modes let you ventilate with enhanced safety, confidence and control.

Thanks to its high-precision ventilator, compact breathing system design and warmer, the Fabius GS premium is an excellent choice for true low flow ventilation. The quick start feature allows manual ventilation in standby mode and during power up for rapid support in emergency situations. Ventilation can be continued even in case of gas supply failure.

COMPREHENSIVE CONTROL

The Fabius GS premium features an integrated, high contrast color monitor which displays all relevant ventilation parameters and curves. The Dräger standard user interface, featuring our intuitive select-adjust-confirm concept makes learning to use the Fabius GS premium very easy. Virtual flow tubes are shown on the high resolution TFT display, and fresh gas consumption data can be recorded and exported for further analysis. An additional mechanical gas flow tube continues to show total gas flow even during a power failure.

CUSTOMIZE YOUR FABIUS

Because every hospital is different, the Fabius GS premium was designed for flexibility. That's why we've included ample mounting space for additional monitors or IT solutions as well as a wide array of useful options. Additional ventilation options and our low flow breathing system kit can extend your ventilation capabilities. Or, add a Dräger Vamos or SCIO monitor for enhanced gas monitoring and control. Standard mounting rails and additional, medical grade power outlet strips let you quickly mount extra Dräger Infinity monitors or even approved third-party equipment. Customize the Fabius GS premium to suit your individual requirements. If your needs change in the future, the Fabius GS premium can change with you. Your Dräger representative will be glad to assist you in choosing the options that are right for you.



ENHANCED ERGONOMICS AND ECONOMICS

The Fabius GS premium was designed around your workflow, not vice versa. The newly designed, highly maneuverable trolley with its convenient central brake makes positioning the Fabius GS premium in the often cramped OR environment easy. Three large drawers provide ample storage space for additional equipment

and consumables, and a large table facilitates comfortable documentation. You have the option to mount the breathing system on either side to better conform to your individual operating room setup. Our DrägerSorb Clic disposable absorbers effectively prevent contact with soda lime and, because you can even exchange them during operation, you can get the most out of every canister.



Advanced information management

The Fabius GS premium is ready for integration into your existing hospital information system. You can then enjoy the benefits of rapid access to the electronic patient file and other findings and even export data to your network – without ever leaving the patient's side.

Additional Dräger Omega or Delta XL monitors can let you maintain a complete, simultaneous overview of ventilation parameters, hemodynamics, diagnostics, laboratory and other parameters. The Fabius GS premium can also let you take advantage of our Infinity monitoring solutions, including our Pick and Go and auto reconfiguration systems for rapid and seamless patient and patient data transfers.



MT-8741-2006

Integrated illumination




D-15974-2009

Total fresh gas flowmeter



MT-2003-2008

CLIC Absorber free
Disposable CO₂ absorber
for Dräger anesthesia
machines

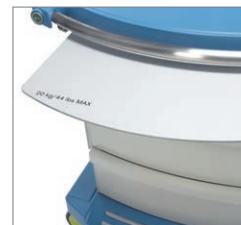

MT-8126-2005

Documentation table

Options and accessoires

A wide array of options and accessories is available to choose from, giving you the power to build exactly the workstation you need:

- Pressure Control, Pressure Support and SIMV/PS ventilation option: Provide ventilation for neonates, pediatrics and adults of practically any acuity level
- Low flow kit (with breathing system warmer and water trap): Avoid condensation build-up in the breathing circuit during low flow ventilation
- External fresh gas outlet
- Mechanical airway pressure gauge
- Dräger Vamos gas monitor
- Auxiliary oxygen flow meter
- Additional RS-232 interface: Additional import and export capability
- Dräger Infinity patient monitors
- Various monitor mounting solutions and a flexible cable management arms
- Additional standard mounting rails with optional rail covers
- Medical grade power outlet strip: Mount and power additional Dräger or third-party equipment, additional monitors and manage device cables
- Pivoting side tray and side basket
- Extendable documentation table: Workflow enhancements that are there when you need them, gone when you don't
- Additional medical-grade lamp: In addition to the built-in LED, an additional lamp can be mounted for extra illumination of the documentation table
- Vaporizer park holder: Holds vaporizers in place when not in use




MT-3850-2008

Power outlet strip
One cable from the power outlet, eliminates tangled cables



MT-5882-2006

Flexible breathing bag extension
for ergonomic positioning of the breathing bag

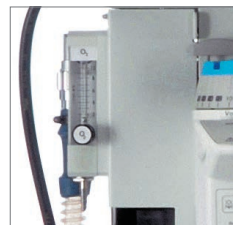

D-354-2010

Auxiliary optional flow meter



MT-3776-2008

Halogen lamp
Additional light source

CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

www.draeger.com

Manufacturer:

Dräger Medical GmbH
Moislinger Allee 53–55
23558 Lübeck, Germany

REGION EUROPE CENTRAL AND EUROPE NORTH

Dräger Medical GmbH
Moislinger Allee 53–55
23558 Lübeck, Germany
Tel +49 451 882 0
Fax +49 451 882 2080
info@draeger.com

REGION EUROPE SOUTH

Dräger Médical S.A.S.
Parc de Haute Technologie d'Antony 2
25, rue Georges Besse
92182 Antony Cedex, France
Tel +33 1 46 11 56 00
Fax +33 1 40 96 97 20
dlmfr-contact@draeger.com

REGION MIDDLE EAST, AFRICA

Dräger Medical GmbH
Branch Office
P.O. Box 505108
Dubai, United Arab Emirates
Tel +971 4 4294 600
Fax +971 4 4294 699
contactuae@draeger.com

REGION ASIA / PACIFIC

Draeger Medical
South East Asia Pte Ltd.
25 International Business Park
#04-27/29 German Centre
Singapore 609916, Singapore
Tel +65 6572 4388
Fax +65 6572 4399
asia.pacific@draeger.com

REGION CENTRAL AND SOUTH AMERICA

Dräger Panama Comercial S. de R.L.
Complejo Business Park,
V tower, 10th floor
Panama City
Tel +507 377 9100
Fax +507 377 9130
contactcsa@draeger.com

CANADA

Draeger Medical Canada Inc.
2425 Skymark Avenue, Unit 1
Mississauga, Ontario, L4W 4Y6
Tel +1 905 212 6600
Toll-free +1 866 343 2273
Fax +1 905 212 6601
Canada.support@draeger.com

As of August 2015:
Dräger Medical GmbH changes
to Drägerwerk AG & Co. KGaA.